


ST70-2/ST60-4 The plate thermo shaker incubator adopts micro-processing techniques combined with PID control, low profile and nice exterior, low noise and easy handling. It is mainly used in mixing and incubating samples in micro-plates, cell culture plates and other standard plates at proper temperature.



Microtiter plates
4 pieces
micro-plates

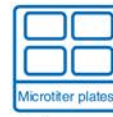

1s-999min
Time Setting



LCD Display
LCD Display

ST70-2



ST60-4



【 Product Features 】

1. LCD display the setting and actual temperature, time and shaking speed simultaneously.
2. Simple touch operation interface, streamlined fuselage, pretty appearance and convenient for cleaning.
3. Both cover and bottom two-sides heating makes sample in each plate hole be heated uniformly.
4. Micro-processor controlling temperature and shaking speed. Good linearity of temperature control. Precise control of shaking speed.
5. Configure 2 plates (for type ST70-2) or 4 plates (for type ST60-4). Beep alarm when operation ends.
6. DC brushless motor drive, long life and free of maintenance.

【 Technical Parameters 】

Model	ST70-2	ST60-4
Speed Range	300 ~1200rpm	300 ~1200rpm
Shaking Orbit	3mm (horizontal)	2mm (horizontal)
Temp. Control Range	R.T.+5℃ ~70℃	R.T.+5℃ ~70℃
Temp. Setting Range	5℃ ~70℃	5℃ ~70℃
Temp. Stability@37℃	±0.2℃	±0.2℃
Platform Temp. Uniformity	≤±0.5℃	≤±0.5℃
Temp. Display Accuracy	0.1℃	0.1℃
Heating Speed Range	≤25min (20℃ to 70℃)	≤25min (20℃ to 70℃)
Time Range	1min ~99h59min or shot mix	1min ~99h59min or shot mix
Sample Capacity	2 Pieces of Plates	4 Pieces of Plates
Voltage	AC 220V / AC 110V, 50/60Hz	AC 220V / AC 110V, 50/60Hz
Power	150W	250W
Fuse	250V, 1A/2A, Φ5x20	250V, 3A/5A, Φ5x20
Dimension (WxDxH)	284 x 264 x 157mm	345 x 310 x 178mm
Net Weight	6.7kgs	10.5kgs



SciLution

บริษัท ไซลิวชั่น จำกัด SciLution Co., Ltd.

Tel : 02 964 8018 | Line ID : @scilution

www.scilution.co.th / www.scilution.net / www.scilution.com